

# NEW ASKIR 230/12V BR

## MAIN APPLICATIONS

Ambulance Car

Emergency Dept.


Home-Care

Patient Transfer

Room / Ward

Tracheostomy

Veterinary

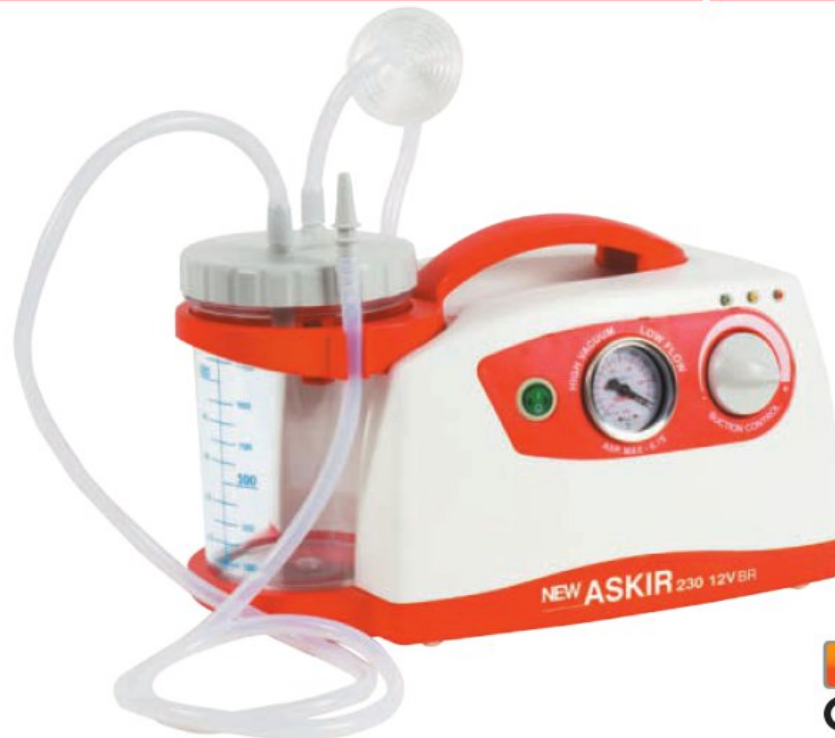

Internal Rechargeable Battery

CE 0123

NEW ASKIR 230/12V BR is an electric suction unit for the aspiration of body liquids, oral, nasal and tracheal aspiration in adults or children. Suitable for emergency applications and for post-operative therapy in home-care. Three different options for operation: AC/DC adapter, rechargeable battery and 12V makes it a versatile suction unit approved for NON-STOP operation without overheating. Provided with visual and acoustic alarm indicating low battery level, the main unit is also equipped with vacuum-meter (bar & kPa) and vacuum regulator. Complies with EN 60601-1-11 "for use in the home healthcare environment" and carries IP21 degree of protection.

## ACCESSORIES INCLUDED

- Liquid Collection Jar with overflow valve system (different options, see below)
- 6x10 Silicone Tubes (autoclavable) + Conical Connector  $\varnothing$  8-9-10mm
- Antibacterial & Hydrophobic Filter (single-patient)
- CH20 Canula
- Power Cord with Schuko plug + 12V Car Adapter

## OPTIONS

- Stand on Castors REF. 27731
  - Carrying bag\* REF. SP 0207
- \*for Unit and Consumables

## TECHNICAL FEATURES

Motor	Oiless and maintenance-free piston pump
Power Feeding	230V-50Hz
	— Internal Rechargeable Battery (12V $\approx$ 2,8 A) - Autonomy 45min - Recharging Time 120÷150min
	— 12V $\approx$ from Car Battery cigarette lighter
ISO 10079-1 Classification	HIGH VACUUM / LOW FLOW
Max Vacuum (adjustable)	-0.75 Bar -75 kPa -563 mmHg
Max free air flow rate	16 l/min
Noise Level	63,0 dB
Duty cycle	Non-stop operation
Weight	Kg 4.90
Size	cm 35 x 21 x 18
Years of Warranty	2
Shipping carton	4
Place of Manufacturing	Italy

Italian  
Medical  
Touch  
MADE IN ITALY

## AVAILABLE CONFIGURATIONS

RE 310200/01	NEW ASKIR 230/12V BR with 1000ml autoclavable jar in polycarbonate
RE 310200/04	NEW ASKIR 230/12V BR with 2000ml autoclavable jar in polycarbonate
RE 310201/01	NEW ASKIR 230/12V BR with 1000ml Flovac <sup>®</sup> disposable liner and reusable container
RE 310200/25	NEW ASKIR 230/12V BR with 2000ml Flovac <sup>®</sup> disposable liner and reusable container
	— Jars in APEC <sup>®</sup> autoclavable up to 143°C available upon request